



Create a clean world by heart



TOP SPEED


Faster than you can imagine

SHENZHEN MARGE TECHNOLOGY SHARING CO.,LTD.
Marge Science and Technology Industrial Park,Shanghenglang Community,
Dalang Street, Longhua District, Shenzhen City,China
Tel: 0755-23219066 23113313 Fax: 0755-23219033
www.margemed.com
Email:tong.chan@margemed.com, sales@margemed.com



Drying object

For all kinds of hard and soft lumen apparatus & instruments, gynecological lumen instruments, ophthalmic precision instruments, dental handpiece, conventional surgical apparatus, instruments and Apparatus not resistant to high temperature

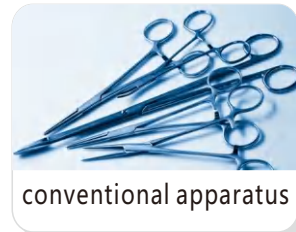

conventional apparatus




ophthalmic instruments



Endoscopic instrument



Orthopedic apparatus

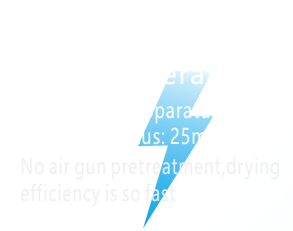

dental apparatus & instrument



medical hose



Anesthetic breathing hose



high temperature resistant rubber hose

Professional solution

Customer needs are our responsibility and mission. In view of the urgent need of large capacity fast low temperature drying equipment is difficult to reach, beauty "Top speed", using a series of innovative technologies, to provide a comprehensive solution for your concerns!



Six advantage



Extreme vacuum technology, so that the boiling point of water is reduced quickly, so that the apparatus can be dried quickly and thoroughly at low temperature



Infrared heating, heating uniformity, fast heating, good drying effect



Humidity detection function, automatically monitor and judge drying effect



Digital display temperature, drying time, vacuum degree



Drying time, reminding function



The whole process overtemperature protection system is safe and reliable

Name: Top speed vacuum Drying cabinet Model :MAG-ZKGZ150

Overall size :775*815*1735mm
Inside groove size 530*695*725mm
Dry shelves: 6
Load :12 baskets
Compatible with basket size:
660*230*50mm
498*230*50mm
480*250*50mm
Drying temperature: normal temperature ~120° C
Cabin body material: stainless steel mirror plate electricity
Power supply:AC380V 3phases 5 wires
Opening mode: double door/single door
Working cycle: 15-25min(optional)



After service system

Afterserive and maintain replacement kits

Top speed vacuum drying
cabinet standard accessories

Top speed vacuum drying cabinet standard after-sales replacement parts				
No	Name	Technology Spec.	Unit	QTY
1	Electric heating plate	220V 500W/pc	pc	2
2	Electric heating tube	220V 300W/pc	pc	2
4	sealing rubber strip	721*551	pc	1
5	PU Tube	Φ 8*Φ6	M	10
6	Fuse Holder	2A	pc	10
9	High temperature hose	Φ45	M	10
10	LED light	3inch white color	pc	1